

# SPECTROSCOUT

PORTABLE ED-XRF SPECTROMETER  
FOR LABORATORY & AT-LINE TESTING



Bringing laboratory-quality  
analysis to the production line

# Fast, reliable, accurate analysis at the point of production

## SPECTROSCOUT

Real-time analytical measurements are essential for monitoring the elemental composition and quality of many processed or manufactured goods.

Users require dependable accuracy and sensitivity at all detection levels, depending on the application, from high-percentage concentrations to trace elements. When results are off-specification, reliable analyses are also key to making timely, data-based process changes.

Many organizations maintain high power wavelength-dispersive XRF or ICP spectrometers in their laboratories to determine elemental compositions. Unfortunately, transporting samples from line to lab wastes precious time. Others locate the spectrometers at the line, risking impaired operation or failure of these sensitive instruments in the production environment.



### SPECTROSCOUT PORTABLE ENERGY DISPERSIVE X-RAY FLUORESCENCE (ED-XRF) ANALYZER

Unlike repackaged handheld X-ray analyzers, SPECTROSCOUT delivers much of the analytical power of top-grade laboratory benchtop instruments but in a durable, compact, easily portable form. SPECTROSCOUT can be calibrated in the lab and then sent straight to the manufacturing floor. It's simple to operate with a little introductory training. And it produces audit-traceable results, right at the point of production – all at surprisingly affordable cost.

#### EXCEPTIONAL PERFORMANCE

SPECTROSCOUT provides outstanding accuracy at concentration levels from ppm to percentages. When measuring high concentrations of several elements, superior sensitivity lets it produce reliable results very quickly. And it gives good light-element performance for powder or solid samples with no need for heavy, costly helium flushing tanks. In fact, the SCOUT is ideal for producers of processed or manufactured goods including chemicals, oil, metal surface coatings, pharmaceuticals, cosmetics, food, feedstock, and more.




#### EXCELLENT EASE-OF-USE

Line personnel find SPECTROSCOUT a snap to operate, without excessive training. Simplified software interacts via a standard computer. Predefined application packages fit the task at hand. Unique iCAL calibration takes just one integrated sample and 5 minutes, reducing recalibration work. And light weight (12 kg/26.46 lb) plus small size (350 x 306 x 306 mm/10.7 x 12.1 x 12.1 in) provide convenient portability and space-saving fit.


# SPECTROSCOUT

Laboratory quality with the right results on the line

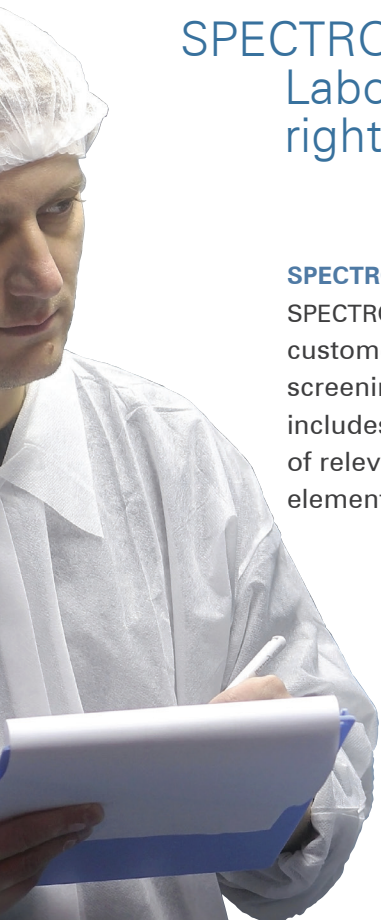

## SPECTROSCOUT COMPACT DESIGN

SPECTROSCOUT offers lab-class advantages with options such as convenient customer-specific calibrations, or advanced unknown-sample extrapolative screening. Superlative standard performance for every SPECTROSCOUT includes good spectral resolution and low limits of detection across a range of relevant materials — from heavy elements such as uranium to light elements such as sodium.

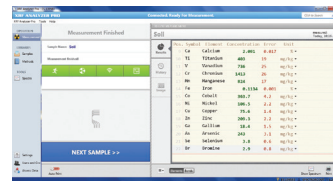

## FAMILIAR INTERFACE

Intuitive design increases ease and productivity without extensive user training.



## PRECISION AND ACCURACY

SPECTROSCOUT provides highly precise measurements from trace to minor to major concentrations of relevant elements. That kind of precision — unmatched by other portable elemental analyzers in its class — enables it to reach performance levels not previously possible.

## COMPETITIVE COST

A SPECTROSCOUT analyzer's purchase price is highly competitive with other portable ED-XRF instruments (while its performance surpasses them). In many cases, its price is also half that of dedicated laboratory instruments. Finally, since its fast, accurate information results in less off-spec product and rework — for greater overall productivity — SPECTROSCOUT promises significant savings for years to come.

## GLOBAL SUPPORT

Industrial operations demand that their instruments be ready for challenging use whenever needed. To ensure that SPECTROSCOUT analyzers are always up to the task, SPECTRO provides the AMECARE Performance Services program.

More than 200 AMECARE service engineers in over 50 countries help ensure optimum performance and extended life for every SPECTROSCOUT instrument. AMECARE's high-value, customized services include proactive maintenance programs; application solutions; access to relevant experts; and instrument-specific training.





# SPECTROSCOUT Specifications



<b>Detector</b>	Silicon drift detector (SDD)
<b>Excitation</b>	X-ray tube, rhodium (Rh) anode, 50 kV max.
<b>DIMENSIONS</b>	
<b>Height</b>	350 mm (13.8 in)
<b>Width</b>	306 mm (12.1 in)
<b>Depth</b>	306 mm (12.1 in)
<b>Weight</b>	12 kg (26.46 lb)
<b>Spectrometer</b>	Integrated PC with touchscreen Bluetooth, wireless LAN, USB interfaces Netbook
<b>Software</b>	SPECTRO XRF Analyzer iCAL (Intelligent Calibration Logic) Data Manager
<b>Accessories</b>	USB cable Consumables
<b>Options</b>	Printer Helium (He) flush and sample spinner Vacuum and sample spinner Battery pack Transport case Video camera

## FULL ANALYZER FAMILY

SPECTRO provides a comprehensive suite of advanced elemental analyzers. ED-XRF analyzers include the high-end SPECTRO XEPOS benchtop analytical instrument, the small spot analyzer SPECTRO MIDEX, the portable yet powerful SPECTROSCOUT and the SPECTRO xSORT handheld XRF spectrometers.

Whatever the product, SPECTRO's more than 30 years of experience in elemental analysis and unparalleled record of technological innovation ensure excellent results.

## www.spectro.com

### GERMANY

SPECTRO Analytical Instruments GmbH  
Boschstrasse 10  
D-47533 Kleve  
Tel. +49.2821.892.0  
spectro.sales@ametek.com

### U.S.A.

SPECTRO Analytical Instruments Inc.  
50 Fordham Rd  
Wilmington 01887, MA  
Tel. +1 800 548 5809  
+1 201 642 3000  
spectro-usa.sales@ametek.com

### CHINA

AMETEK Commercial  
Enterprise (Shanghai) CO., LTD.  
Part A1, A4 2nd Floor Building No. 1 Plot Section  
No. 526 Fute 3rd Road East; Pilot Free Trade Zone  
200131 Shanghai  
Tel. +86.400.022.7699  
spectro-china.sales@ametek.com



### Subsidiaries:

- **FRANCE:** Tel. +33.1.3068.8970, spectro-france.sales@ametek.com
- **GREAT BRITAIN:** Tel. +44.1162.462.950, spectro-uk.sales@ametek.com
- **INDIA:** Tel. +91.22.6196.8200, sales.spectroindia@ametek.com
- **ITALY:** Tel. +39.02.94693.1, spectro-italy.sales@ametek.com
- **JAPAN:** Tel. +81.3.6809.2405, spectro-japan.info@ametek.co.jp
- **SOUTH AFRICA:** Tel. +27.11.979.4241, spectro-za.sales@ametek.com

SPECTRO operates worldwide and is present in more than 50 countries. For SPECTRO near you, please visit [www.spectro.com/worldwide](http://www.spectro.com/worldwide)  
© 2015 AMETEK Inc., all rights reserved, subject to technical modifications • F-15, Rev.1 • 80902730 • Photos: SPECTRO, Adobe Stock  
Registered trademarks of SPECTRO Analytical Instruments GmbH • **SPECTRO**: USA (3,645,267); EU (005673694); "SPECTRO": EU (009693763);  
"SPECTROSCOUT": EU (10705424); USA (4,394,468); Japan (1141441); South Korea (1141441); China (1141441)